


Superior Irrigation and Evacuation for Endoscopy




Clear Upper GI Bleeding

Fast, convenient, reliable, clear view
Seeing means completing

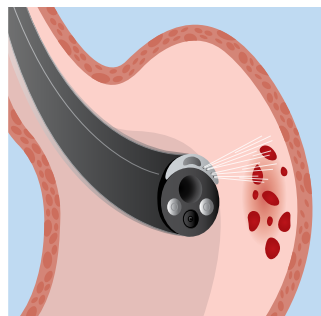

Concurrent Working Channel:
ClearPath Upper GI handles all irrigation and evacuation. Endoscope working channel always available.

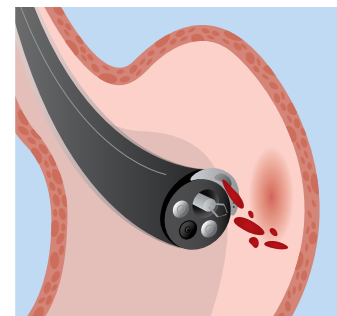

Evacuate Larger Clots:
Ordinary endoscope evacuation cross section 8-12 mm².
ClearPath Upper GI evacuation cross section 27 mm²

Superior Irrigation and Evacuation for Endoscopy



Clear Upper GI Bleeding

Concurrent working channel

Eliminates use of endoscope for evacuation.
Endoscope working channel always available.

Evacuate larger clots

Ordinary endoscope evacuation cross section 8-12 mm².
ClearPath evacuation cross section 27 mm².

Rapid deployment

From open kit to procedure-ready.

Fits all common upper GI endoscopes

Head attaches over distal end. Irrigation/evacuation tubes attach along length. Operates by standard foot pedals.

Product Specification	
Irrigation pressure	1.5 - 2.5 bar
Evacuation rate	1800 ml/min
Evacuation cross section	27 mm ²
Use with	ClearPath Controller

EasyGlide Ltd.






30 Ha'Ella St., Kfar Truman 73150, Israel Fax/Phone: 972.77.6201003
info@easy-glide.com • www.easy-glide.com

ClearPath is a trademark of Easy-Glide Ltd.

©2011 Easy-Glide, Ltd. All rights reserved

CE 0470